

MXR Podoblock™

Innovative vet  solutions
through experience.





WE DON'T JUST SELL GOOD PRODUCTS
THEY'RE **WORKS OF ART**, REALLY



Phil wants
to order
12 Jixtrs.
Call back!



Nico Gorter
(+31) 6 53 79 42 82

Alex Berends
(+31) 6 51 82 81 17



Equine Gazette

MXR's Stat-X Espléndido comes with all-new 3D Head

Præsent vestibulum molestie lacus. Aenean nonummy hendrerit mauris. Phasellus porta. Fusce suscipit Donec eget tellus non erat lacinia fermentum. Donec in velit vel ipsum auctor pulvinar. Proin ullamcorper urna et felis. Fusce suscipit ridiculus mus. Nulla dui. Fusce feugiat malesuada odio. Morbi nunc Duis ultricies pharetra magna.



0031 50 8200 257

MXR Podoblock B.V.

Nico Gorter / Alex Berends

Company address:
Esweg 31
9475 TM Midlaren
The Netherlands

Tel: Nico 0031-653 794 282
Alex 0031-651 828 117
Office 0031-508 200 257
Fax 0031-508 200 212

Web: www.podoblock.com
e-Mail: info@mxr.nl
e-Mail: sales@mxr.nl

Chamber of commerce: 58189998
VATnr.: NL852918227B01



Welcome to our world!

Hereby we proudly present the new product catalogue from MXR Podoblock. We are pleased to offer you a look into our complete product line. All our products have been invented in order for you to make quick and safe, high-quality X-rays. Every piece of equipment has been innovated, designed, built and tested by ourselves. They are purely based on experience, earned during 15 years of working daily with equine X-rays.



As of now most of our products will be fitted with a NFC-tag. This 'Near Field Communication'-tag will give you and your mobile phone the opportunity to receive additional practical information on that specific product, like instructional videos or safety tips and tricks for optimal working situations. NFC works similarly to the old QR-codes, just simply hold your mobile phone close to the tag and all the information needed will appear on your screen.

Of course we remain active on , ,  and , so follow us or link up and you will be the first in line to get the latest news and product updates. We work with a worldwide and direct network to stay in contact with veterinarians and distributors. They are our source for the best ideas and new problems to solve. This way we can always offer you the latest information and come in first with new products. Please challenge us with your ideas or problems and we may very well come up with a solution.







All products you find inside are made from the highest quality materials and are manufactured in Holland. This is of course for daily use and to last for years. MXR Podoblock supplies products for every X-ray system, both for standard and, if desired, special tailor-made products. It doesn't matter if you use a conventional, analogue, digital CR or the latest digital DR flat-panel system, we offer solutions for all of them.

We supply our products worldwide through a network of distributors or, if needed, directly to you. Feel free to roam around at our website for a total list of distributors. They supply all of our products. As you can see in this catalogue, we have standard solutions for almost everything in our product line. Our latest catalogue will always contain novelties and innovations we are working on. If you are not sure about something or you can't find the product you are looking for in the catalogue, please contact us directly and we will surely offer you the solution you need!

We hope to welcome you into our world,

Nico & Alex





| Pictures | Product Description | Specifications | Article number |
|--|---|---|--|
|  | <p>Podoblock.</p> <p>Aluminium. Incl. 1 Topblock. 4 Rubber Feet. LM / AP / 45° Hands-free system.</p> | <p>In combination with our Podoholder, Podo-Prot.box or Cowboybox this unique system guarantees maximum radiation distance and no fingers on the X-rays.</p> <p>The height from ground to tip off hoof in Podoblock: 90 mm 120 mm 170 mm</p> | <p>PB 01AL090 PB 01AL120 PB 01AL170</p> |
|  | <p>Bottom plate.</p> <p>Synthetic PE500. White.</p> | <p>Option to create a little sliding effect. Allows for small adjustments even with pressure of the horse's feet placed on the Podoblock. Thickness of the plate is 10mm.</p> | <p>PV 01PE500</p> |
|  | <p>Topblock.</p> <p>Synthetic PE500. Black. Angle diversity from 45 ° to 60 °.</p> <p>Topblock Wedge . To be placed on our standard Topblock.</p> <p>Topblock Flatplate.</p> <p>Make LM / AP / 45 ° with the foot standing flat or upright pedal</p> | <p>The Topblock is easily interchangeable.</p> <p>45 ° 50 ° 55 ° 60 ° 60 ° Upright pedal.</p> <p>10° extra angle.</p> <p>Flatplate AP Flatplate LM & AP Flatplate LM & AP & 45 °</p> | <p>TB 0250045 TB 0250050 TB 0250055 TB 0250060 TBD 0250060</p> <p>TBW 0250010</p> <p>TBF 02500AP TBF 02500LMAP TBF 0250045</p> |

| | | | |
|--|---|---|---|
|   | <p>Podoholder.</p>  <p>To be used in combination with Podoblock and Combibox. Protects the cassette. Secures always the right positions for standard hoof exams.</p> <p>Podoholder Lead. Back-side with lead integrated against X-ray scatter.</p> <p>Podoholder Grid. Can be provided to place your own or our Grid in the compartment. This way the Grid is always protected and exactly in the centre.</p> | <p>18 x 24 18 x 24 Horizontal 24 x 30 24 x 30 Horizontal 35 x 43 Horizontal 8" x 10" 8" x 10" Horizontal 10" x 12" 10" x 12" Horizontal 10" x 14" 10" x 14" Horizontal</p> <p>18 x 24 Lead 24 x 30 Lead 8" x 10" Lead 10" x 12" Lead 10" x 14" Lead</p> <p>18 x 24 Grid 24 x 30 Grid 8" x 10" Grid 10" x 12" Grid 10" x 14" Grid</p> | <p>PH 03182418 PHB 03182418 PH 03243018 PHB 03243018 PHB 03354318 PH 03232818 PHB 03232818 PH 03283318 PHB 03283318 PH 03101418 PHB 03101418</p> <p>PL 03182418 PL 03243018 PL 03232818 PL 03283318 PL 03101418</p> <p>PR 03182418 PR 03243018 PR 03232818 PR 03283318 PR 03101418</p> |
|   | <p>Distance Cassette Holder.</p> <p>Invented and built to create maximum distance from the radiation as well from the horse. 3D-adjustable telescopic stick with only one knob to secure. Super strong durable lightweight Polycarbonate.</p>  <p>Distance Cassette Holder with foot option.</p> <p>Same as our standard Dist.cass.holder but with a mountable and in height adjustable foot.</p> <p>Available in 2 lengths: 300 mm for Hock. 1200mm for Back / Neck. Please recommend by ordering the desired length of the foot.</p> | <p>18 x 24 18 x 24 Fetlock 24 x 30 35 x 43 35 x 43 Horizontal 8" x 10" 8" x 10" Fetlock 10" x 12" 10" x 14" 20 x 30 20 x 40 30 x 40</p> <p>18 x 24 Stative 24 x 30 Stative 35 x 43 Stative 35 x 43 Stative, Horizontal 8" x 10" Stative 10" x 12" Stative 10" x 14" Stative 20 x 30 Stative 20 x 40 Stative 30 x 40 Stative</p> | <p>AC 04182418 ACK 04182418 AC 04243018 AC 04354318 ACB 04354318 AC 04232818 ACK 04232818 AC 04283318 AC 04101418 AC 04203018 AC 04204018 AC 04304018</p> <p>ACS 04182418 ACS 04243018 ACS 04354318 ACBS 04354318 ACS 04232818 ACS 04283318 ACS 04101428 ACS 04203018 ACS 04204018 ACS 04304018</p> |










Note: These items above are only suitable for use with Analog and CR cassettes.





| | | | |
|---|---|---|--|
|  | <p>Skylinebox. </p> <p>Most effective product on the market to take high quality Skyline X-rays.</p> <p>Makes taking Skyline exams simple because of the special built-in angle.</p> <p>After the X-ray exposure possible to reload the cassette without replacing the horse.</p> | <p>Available for all sizes Analogue and CR cassettes. Rubber feet. Polycarbonate 10mm 3M Antislip. Skyline Angle 5°.</p> <p>18 x 24 24 x 30 8" x 10" 10" x 12" 10" x 14"</p> | <p>SK 05182420 SK 05243020 SK 05232820 SK 05283320 SK 05101420</p> |
| <p>Combibox turned to Skyline</p>  <p>Combibox with Podoholder for LM / AP</p>  <p>Combibox LM / AP & Dummy</p>  | <p>Combibox. </p> <p>Above all the features of the Skylinebox when turned over, also to be used for LM / AP flat and upright pedal.</p> <p>With standard hands-free system with the optional Podoholder. 1 box that features making Skyline, LM & AP X-rays of the foot.</p> <p>For the optimal solution use it in combination with the Dummy box to secure that both feet are standing on the same height and the horse is standing comfortably. Makes it possible to take X-rays from both feet without moving the horse.</p> | <p>Available for all sizes cassette, Analogue and CR cassettes. Rubber feet. Polycarbonate 10 & 18mm Height 70 mm. 3M Antislip. Skyline Angle 5°. Hands free system with Podoholder: LM or LM & AP</p> <p>18 x 24 LM 18 x 24 LM / AP 18 x 24 LM / AP & Dummy</p> <p>24 x 30 LM 24 x 30 LM / AP 24 x 30 LM / AP & Dummy</p> <p>8" x 10" LM 8" x 10" LM / AP 8" x 10" LM / AP & Dummy</p> <p>10" x 12" LM 10" x 12" LM / AP 10" x 12" LM / AP & Dummy</p> <p>10" x 14" LM 10" x 14" LM / AP 10" x 14" LM / AP & Dummy</p> <p>Dummy LM</p> | <p>CO 061824081 CO 061824082 CO 061824083</p> <p>CO 062430081 CO 062430082 CO 062430083</p> <p>CO 062328081 CO 062328082 CO 062328083</p> <p>CO 062833081 CO 062833082 CO 062833083</p> <p>CO 061014081 CO 061014082 CO 061014083</p> <p>CO 06DUM081</p> |

Note: These items above are only suitable for use with Analog and CR cassettes.

| | | | |
|---|--|--|--|
|  | <p>Knee Holder.</p> <p>This new product is a simple to use distance holder that create a maximum distance and safety.</p> <p>Is standard supplied with a 400mm stick that can be adjusted to the left or right in a split second.</p> <p>The material for the Knee Holder used is Polycarbonate 3mm.</p> |  <p>For all Analogue and CR cassettes.</p> <p>18 x 24 24 x 30 8" x 10" 10" x 12" 10" x 14"</p> | <p>KH 04182418 KH 04243018 KH 04232818 KH 04283318 KH 04101418</p> |
|  | <p>Dental Protection box.</p> <p>Protects your cassette to safely take intraoral X-rays of the incisors and canines and LM and AP of the head.</p> <p>The material for the Protectionbox used is Polycarbonate 6mm.</p> <p>The mountable 400mm long stick, can be turned 360 °.</p> <p>This feature guarantees a minimum of radiation for the person holding the cassette.</p> <p>Easy to use with horses because you can fully concentrate on the position of the cassette.</p> <p>The Dental protection box can also be used together with our Denfinity.</p> |  <p>For all Analogue and CR cassettes.</p> <p>18 x 24 24 x 30 8" x 10" 10" x 12" 10" x 14"</p> | <p>MB 04182416 MB 04243016 MB 04232816 MB 04283316 MB 04101416</p> |

Note: These items above are only suitable for use with Analog and CR cassettes.

| | | | |
|---|--|---|--|
|     | <p>Denfinity.</p> <p>This new light weight portable and collapsible stainless steel stand offers you the freedom to make high quality Skull X-rays.</p> <p>As a start you use the Denfinity stand. This fully stainless steel base plate and vertical telescopic pole guarantees always a secure positioning and right height adjustment.</p> <p>The Denfinity stand allows you to place and replace your required option in a split seconds. You can choose between the Aluminium 3D swivel-tilt joint in combination with our Dentalbox or Protectionbox.</p> <p>The LM-Head support and the fully radiolucent Headrest can be placed directly on the Denfinity.</p> <p>The Denfinity allows you to make all Skull and Jaw X-rays hands free. There are no extra persons needed to take high quality images and create a save working environment.</p> <p>In progress to order the optional Protectionbox, Podoholder or Dentalbox place with your order the required cassette measurements you use. If you have any question about this, contact your local supplier or MXR.</p> | <p>Denfinity. Base plate: 350x405x90 mm Vertical pole: 960 mm Vertical travel: 600 mm Material: Stainless steel</p>  <p>3D-Joint. Aluminium 3D swivel-tilt joint Used in combination with Dentalbox or Protectionbox.</p>  <p>LM-Head support. Material: Polycarbonate 10mm Swivel range: 180° Universal adapter for fixation Podoholder or Protectionbox</p>  <p>Headrest. Material: Polycarbonate 10mm Full radiolucent, no metal parts in detection area at all.</p>  <p>Dentalbox. Material: Polycarbonate 6mm Delivered with 400mm stick for extra user options.</p>  <p>Protectionbox. Material: Polycarbonate 3mm Fully closed cassette or DR-panel Protectionbox.</p> | <p>DSJ 1050401</p> <p>DSJ 10DNL3D</p> <p>DSJ 10LMHEAD</p> <p>DSJ 10HREST</p> <p>DSJ 10DENT</p> <p>DSJ 10PROT</p> |
|---|--|---|--|

| | | | |
|--|---|---|----------------------------------|
|  | <p>The "Juliuster"</p>  <p>This new product is designed, tested and approved in collaboration with European Equine Dental Consultancy EEDC by Dr. Julius W.E. Peters.</p> <p>Easy to use device to offset the lower jaw to the left or right for open mouth views with offset technique.</p> <p>Great to investigate the roots and the occlusion surface of the cheek teeth.</p> <p>The Dorso-Ventral views can be made with 5 cm offset to the left and right.</p> <p>Delivered with full Biothane strap and metal buckle that goes from Shet to Shire.</p> | <p>Material: HMPE Offset: 5 cm left / right Adjusted in 12mm steps Strap: 148cm Biothane</p> | <p>DRS 1011560</p> |
|  | <p>Dental X-Block.</p> <p>Can be used directly between the incisors to create quickly 60, 80 or 100 mm space.</p> | <p>Material: Rubber 95°Shore Size: 60 x 80 x 100 mm.</p> | <p>DXB 106810</p> |
|  | <p>Sørender.</p> <p>Be used directly between the scissor teeth in combination with full anaesthetic with spare Ø30 or Ø50 mm hole for tube transit. Standard available in Medium and Large.</p> | <p>Material: Polypropylene Strap from Biothane Medium: 100x100x50 / Ø30 Large: 120x115x55 / Ø50</p> | <p>AB 1010030 AB 1012050</p> |

DR Flat panel solutions.

MXR Podoblock has a standard solution for every DR flat panel. It doesn't matter if you use the 'newest' or the 'oldest' model on the market. Because the DR flat panels are expensive and more delicate than a CR cassette, you definitely need an good protection for the panel itself.











We provide a complete and professional range of products to accomplish this. The basic thing you need is a Protectionbox or Sleeve, our Cowboy box "all in one" is the optimal solution.

Of course we provide DR Skyline tunnels to make navicular Skyline exams of the hoof. You can also upgrade from Skylinebox to the optimal solution: our DR-Combibox.

Please ask your local supplier of the DR flat panel about our products.

For a complete list of all DR flat panel manufacturers and suppliers please see the next pages.

According to this list you can decide what MXR Podoblock article number needs to be ordered.
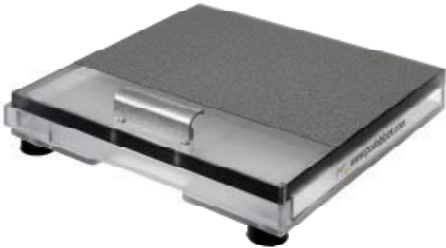
| | |
|---|---|
|  | <p>DR-Protectionbox. </p> <p>Light-weight protection box. Handheld.</p> <p>Optional: - Grid Compartment (tailor-made) - Greppstr 2.1</p> |
|  | <p>DR-Podo-Protectionbox. </p> <p>Protectionbox with 5mm adaptors to place in our Podoblock and Combibox. Guarantees the same direction and angle with hoof examinations at all times.</p> <p>Optional: - Grid Compartment (tailor-made) - Greppstr 2.1</p> |
|  | <p>DR-Cowboybox. </p> <p>Protectionbox, Podo-Protectionbox and Remote cassette holder all in one.</p> <p>Standard with a mountable telescopic grip which can be turned in 3D positions to get extra distance.</p> <p>Mountable and height adjustable support foot 300mm</p> <p>Optional: - Grid Compartment (tailor-made) - Greppstr 2.1</p> |
|  | <p>Greppstr 2.1 </p> <p>This new universal grip handle can be placed in 2 seconds on almost every Protectionbox, Podo-Protectionbox and Cowboybox and also on several DR-Sleeve's.</p> <p>This way you can choose to place your Greppstr 2.1 where you want or leave it away if you want to use a plane box.</p> <p>Note: If you already have a product from us and want to add our Greppstr to it, please contact us and we can arrange it.</p> |
|  | <p>DR-Skylinebox. </p> <p>Ideal box for standard Skyline exams of the foot.</p> <p>Made of 20mm Polycarbonate.</p> <p>Because of the special angle the horse is standing comfortably and there is more space for the X-ray generator.</p> <p>The top side is supplied with a radiolucent anti-slip surface and the door is secured by a lock.</p> |

DR-Combibox LM with optional Podo-Protectionbox



DR-Combibox LM / AP with optional Podo-Protectionbox



DR-Combibox LM.

The Skyline side of the Combibox has the same features as the original DR-Skylinebox.

Because of the special angle the horse is standing comfortably and there is more space for the X-ray generator.

If the DR-Skyline box is turned over it has the extra features to make easy LM and AP flat views of the foot.

The special adapters on the DR-Combibox allows you to take "hands-free" X-rays in combination with our Podo-Prot.box or DR-Cowboybox.

This way the DR flat panel is protected and always in the right position.



DR-Combibox LM / AP.

There is the possibility to choose between only the LM adapter or, your best option of course, LM and AP adapters.

This way it is not necessary to replace the foot for both exams.

DR-Combibox SET.

Both standing feet on the same height is advised with LM and AP examinations. That's why our Combibox Dummy can also be supplied with the DR-Combibox.

The Combibox Set is included:

DR-Combibox LM / AP

Combibox Dummy LM



DR-Cowboy Package.

The optimal combined package of DR-Products to make all X-rays needed safely and in a standardized fashion.

This is what you need, for mobile or stationary use in the clinic.

The package deal includes:

Podoblock

DR-Cowboybox

DR-Combibox LM / AP

Combibox Dummy LM



DR-Sleeve.









Soft cover Sleeve made of PVC with Velcro closing.

Front side is full radiolucent with marked detection area.

Inside foam at the edges and corners for extra protection.

Depending on the DR system you use the Sleeve comes with a mounted grip handle or with our [Grepstr 2.1](#)

Note: The soft cover DR-Sleeve is not suitable to protect your panel against impact force !
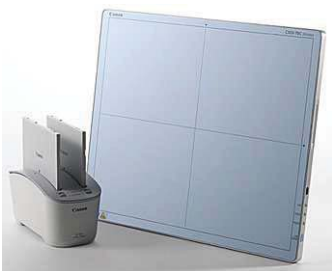
| 1. DR flat panel | 2. Type / Name | 3. Product type | 4. Article number |
|--|--|--|---|
| Canon  | Canon CXDI-31 | DR-Protectionbox DR-Skylinebox DR-Combibox LM DR-Combibox LM / AP DR-Combibox SET DR-Podo-Protectionbox DR-Cowboybox DR-Cowboy Package DR-Sleeve | DR31 03323330 DR31 05323330 DR31 06323330 1 DR31 06323330 2 DR31 06323330 S DR31 07323330 DR31 08323330 DR31 08323330 P DR31 10323330 |
| Canon  | Canon CXDI-50C Canon CDXI-50G Scope-XL XDR1 XDR2 | DR-Protectionbox DR-Skylinebox DR-Combibox LM DR-Combibox LM / AP DR-Combibox SET DR-Podo-Protectionbox DR-Cowboybox DR-Cowboy Package DR-Sleeve | DR50 03514528 DR50 05514528 DR50 06514528 1 DR50 06514528 2 DR50 06514528 S DR50 07514528 DR50 08514528 DR50 08514528 P DR50 10514528 |
| Canon  | Canon CXDI-55C Canon CXDI-55G Sedecal EPX-2 | DR-Protectionbox DR-Skylinebox DR-Combibox LM DR-Combibox LM / AP DR-Combibox SET DR-Podo-Protectionbox DR-Cowboybox DR-Cowboy Package DR-Sleeve | DR55 03484815 DR55 05484815 DR55 06484815 1 DR55 06484815 2 DR55 06484815 S DR55 07484815 DR55 08484815 DR55 08484815 P DR55 10484815 |
| Canon  | Canon CXDI-60C Canon CXDI-60G Eklin Mark IIG Eklin Mark IIIG Scope-XS XDR2 Sedecal EPX-1 | DR-Protectionbox DR-Skylinebox DR-Combibox LM DR-Combibox LM / AP DR-Combibox SET DR-Podo-Protectionbox DR-Cowboybox DR-Cowboy Package DR-Sleeve | DR60 03343824 DR60 05343824 DR60 06343824 1 DR60 06343824 2 DR60 06343824 S DR60 07343824 DR60 08343824 DR60 08343824 P DR60 10343824 |
| Canon  | Canon CXDI-70C Canon CXDI-701C Canon CXDI-701G Sedecal WEPX-2 | DR-Protectionbox DR-Skylinebox DR-Combibox LM DR-Combibox LM / AP DR-Combibox SET DR-Podo-Protectionbox DR-Cowboybox DR-Cowboy Package DR-Sleeve | DR70 03384617 DR70 05384617 DR70 06384617 1 DR70 06384617 2 DR70 06384617 S DR70 07384617 DR70 08384617 DR70 08384617 P DR70 10384617 |
| Canon  | Canon CXDI-80C Canon CXDI-801C Canon CXDI-801G Eklin Mark 1114CW Sedecal WEPX-1 | DR-Protectionbox DR-Skylinebox DR-Combibox LM DR-Combibox LM / AP DR-Combibox SET DR-Podo-Protectionbox DR-Cowboybox DR-Cowboy Package DR-Sleeve | DR80 03383015 DR80 05383015 DR80 06383015 1 DR80 06383015 2 DR80 06383015 S DR80 07383015 DR80 08383015 DR80 08383015 P DR80 10383015 |


| | | | |
|--|---|---|---|
| <p>Sound-Eklin</p>  | <p>Eklin Mark II Idaho Eklin Mark III Idaho</p> | <p>DR-Protectionbox DR-Skylinebox DR-Combibox LM DR-Combibox LM / AP DR-Combibox SET DR-Podo-Protectionbox DR-Cowboybox DR-Cowboy Package DR-Sleeve</p> | <p>DREK 03IDAHO DREK 05IDAHO DREK 06IDAHO 1 DREK 06IDAHO 2 DREK 06IDAHO S DREK 07IDAHO DREK 08IDAHO DREK 08IDAHO P DREK 10IDAHO</p> |
| <p>Sound-Eklin</p>  | <p>Eklin Mark II Texas Eklin Mark III Texas</p> | <p>DR-Protectionbox DR-Skylinebox DR-Combibox LM DR-Combibox LM / AP DR-Combibox SET DR-Podo-Protectionbox DR-Cowboybox DR-Cowboy Package DR-Sleeve</p> | <p>DREK 03TEXAS DREK 05TEXAS DREK 06TEXAS 1 DREK 06TEXAS 2 DREK 06TEXAS S DREK 07TEXAS DREK 08TEXAS DREK 08TEXAS P DREK 10TEXAS</p> |
| <p>Vatech</p>  | <p>Vatech 1210P Samsung 1210P Leonardo DR 1210Pvet Cuattro Uno EQ Slate 1012 Cuattro Uno EQ 1012 XDR2-K, XDR2-M, XDR2-S meX+ 1210P Sedecal EPX-V Eklin Mark 1210C</p> | <p>DR-Protectionbox DR-Skylinebox DR-Combibox LM DR-Combibox LM / AP DR-Combibox SET DR-Podo-Protectionbox DR-Cowboybox DR-Cowboy Package DR-Sleeve</p> | <p>DRSX 03514528 DRSX 05514528 DRSX 06514528 1 DRSX 06514528 2 DRSX 06514528 S DRSX 07514528 DRSX 08514528 DRSX 08514528 P DRSX 10514528</p> |
| <p>Vatech</p>  | <p>1012WCA Wireless Eklin Mark 1210CW Sprint™ Ultralight DR XDR3 meX+ 1012WCA</p> | <p>DR-Protectionbox DR-Skylinebox DR-Combibox LM DR-Combibox LM / AP DR-Combibox SET DR-Podo-Protectionbox DR-Cowboybox DR-Cowboy Package DR-Sleeve</p> | <p>DRSX 031012WCA DRSX 051012WCA DRSX 061012WCA 1 DRSX 061012WCA 2 DRSX 061012WCA S DRSX 071012WCA DRSX 081012WCA DRSX 081012WCA P DRSX 101012WCA</p> |
| <p>Vatech</p>  | <p>1417PCA / PGA 1417WCA / WGA Wireless</p> | <p>DR-Protectionbox DR-Skylinebox DR-Combibox LM DR-Combibox LM / AP DR-Combibox SET DR-Podo-Protectionbox DR-Cowboybox DR-Cowboy Package DR-Sleeve</p> | <p>DRSX 031417PGA DRSX 051417PGA DRSX 061417PGA 1 DRSX 061417PGA 2 DRSX 061417PGA S DRSX 071417PGA DRSX 081417PGA DRSX 081417PGA P DRSX 101417PGA</p> |
| <p>Vatech</p>  | <p>1417WCC/WGC Wireless</p> | <p>DR-Protectionbox DR-Skylinebox DR-Combibox LM DR-Combibox LM / AP DR-Combibox SET DR-Podo-Protectionbox DR-Cowboybox DR-Cowboy Package DR-Sleeve</p> | <p>DRSX 031417WGC DRSX 051417WGC DRSX 061417WGC 1 DRSX 061417WGC 2 DRSX 061417WGC S DRSX 071417WGC DRSX 081417WGC DRSX 081417WGC P DRSX 101417WGC</p> |

| | | | |
|---|---|---|--|
| <p>Carestream</p>  | <p>DRX-1 Wireless Examion Flexible FPS</p> | <p>DR-Protectionbox DR-Skylinebox DR-Combibox LM DR-Combibox LM / AP DR-Combibox SET DR-Podo-Protectionbox DR-Cowboybox DR-Cowboy Package DR-Sleeve</p> | <p>DRKO 03384617 DRKO 05384617 DRKO 06384617 1 DRKO 06384617 2 DRKO 06384617 S DRKO 07384617 DRKO 08384617 DRKO 08384617 P DRKO 10384617</p> |
| <p>Varian 2520 0° connector</p>  | <p>Paxscan 2520E+ Leonardo DR2520+ Protec Rapixx 2520V</p> | <p>DR-Protectionbox DR-Skylinebox DR-Combibox LM DR-Combibox LM / AP DR-Combibox SET DR-Podo-Protectionbox DR-Cowboybox DR-Cowboy Package DR-Sleeve</p> | <p>DRLE 03234528 DRLE 05234528 DRLE 06234528 1 DRLE 06234528 2 DRLE 06234528 S DRLE 07234528 DRLE 08234528 DRLE 08234528 P DRLE 10234528</p> |
| <p>Varian 2520 90° connector</p>  | <p>Varian Paxscan 2520E+ Cuattro Uno EQ Sound Eclin truDR Protec Rapixx 2520V</p> | <p>DR-Protectionbox DR-Skylinebox DR-Combibox LM DR-Combibox LM / AP DR-Combibox SET DR-Podo-Protectionbox DR-Cowboybox DR-Cowboy Package DR-Sleeve</p> | <p>DRCU 03234528 DRCU 05234528 DRCU 06234528 1 DRCU 06234528 2 DRCU 06234528 S DRCU 07234528 DRCU 08234528 DRCU 08234528 P DRCU 10234528</p> |
| <p>Fuji FDR</p>  | <p>D-EVO C24i Wireless</p> | <p>DR-Protectionbox DR-Skylinebox DR-Combibox LM DR-Combibox LM / AP DR-Combibox SET DR-Podo-Protectionbox DR-Cowboybox DR-Cowboy Package DR-Sleeve</p> | <p>DRFU 03C24i DRFU 05C24i DRFU 06C24i 1 DRFU 06C24i 2 DRFU 06C24i S DRFU 07C24i DRFU 08C24i DRFU 08C24i P DRFU 10C24i</p> |
| <p>Fuji FDR</p>  | <p>D-EVO G35</p> | <p>DR-Protectionbox DR-Skylinebox DR-Combibox LM DR-Combibox LM / AP DR-Combibox SET DR-Podo-Protectionbox DR-Cowboybox DR-Cowboy Package DR-Sleeve</p> | <p>DRFU 03DEVO DRFU 05DEVO DRFU 06DEVO 1 DRFU 06DEVO 2 DRFU 06DEVO S DRFU 07DEVO DRFU 08DEVO DRFU 08DEVO P DRFU 10DEVO</p> |
| <p>Fuji FDR</p>  | <p>D-EVO C35i Wireless D-EVO II C35 Wireless</p> | <p>DR-Protectionbox DR-Skylinebox DR-Combibox LM DR-Combibox LM / AP DR-Combibox SET DR-Podo-Protectionbox DR-Cowboybox DR-Cowboy Package DR-Sleeve</p> | <p>DRFU 03DEVOW DRFU 05DEVOW DRFU 06DEVOW 1 DRFU 06DEVOW 2 DRFU 06DEVOW S DRFU 07DEVOW DRFU 08DEVOW DRFU 08DEVOW P DRFU 10DEVOW</p> |

| | | | |
|--|---|---|--|
| <p>Idexx</p>  | <p>Examion DR810 Idexx Equiview Idexx DR 810 FDXD 810</p> | <p>DR-Protectionbox DR-Skylinebox DR-Combibox LM DR-Combibox LM / AP DR-Combibox SET DR-Podo-Protectionbox DR-Cowboybox DR-Cowboy Package DR-Sleeve</p> | <p>DRID 03334225 DRID 05334225 DRID 06334225 1 DRID 06334225 2 DRID 06334225 S DRID 07334225 DRID 08334225 DRID 08334225 P DRID 10334225</p> |
| <p>DR Tech</p>  | <p>Flaatz 560</p> | <p>DR-Protectionbox DR-Skylinebox DR-Combibox LM DR-Combibox LM / AP DR-Combibox SET DR-Podo-Protectionbox DR-Cowboybox DR-Cowboy Package DR-Sleeve</p> | <p>DRFL 03384617 DRFL 05384617 DRFL 06384617 1 DRFL 06384617 2 DRFL 06384617 S DRFL 07384617 DRFL 08384617 DRFL 08384617 P DRFL 10384617</p> |
| <p>Thales</p>  | <p>Pixium Portable 2430EZ Leonardo DR Mini</p> | <p>DR-Protectionbox DR-Skylinebox DR-Combibox LM DR-Combibox LM / AP DR-Combibox SET DR-Podo-Protectionbox DR-Cowboybox DR-Cowboy Package DR-Sleeve</p> | <p>DRPI 032430EZ DRPI 052430EZ DRPI 062430EZ 1 DRPI 062430EZ 2 DRPI 062430EZ S DRPI 072430EZ DRPI 082430EZ DRPI 082430EZ P DRPI 102430EZ</p> |
| <p>PerkinElmer</p>  | <p>XRpad™ 4336</p> | <p>DR-Protectionbox DR-Skylinebox DR-Combibox LM DR-Combibox LM / AP DR-Combibox SET DR-Podo-Protectionbox DR-Cowboybox DR-Cowboy Package DR-Sleeve</p> | <p>DRPE 03XR4336 DRPE 05XR4336 DRPE 06XR4336 1 DRPE 06XR4336 2 DRPE 06XR4336 S DRPE 07XR4336 DRPE 08XR4336 DRPE 08XR4336 P DRPE 10XR4336</p> |
| <p>Konica Minolta</p>  | <p>AeroDR 1417S AeroDR 1417HQ Wireless</p> | <p>DR-Protectionbox DR-Skylinebox DR-Combibox LM DR-Combibox LM / AP DR-Combibox SET DR-Podo-Protectionbox DR-Cowboybox DR-Cowboy Package DR-Sleeve</p> | <p>DRKM 03AERO DRKM 05AERO DRKM 06AERO 1 DRKM 06AERO 2 DRKM 06AERO S DRKM 07AERO DRKM 08AERO DRKM 08AERO P DRKM 10AERO</p> |
| <p>Konica Minolta</p>  <p>1012HQ 1417HQ 1717HQ</p> | <p>AeroDR 1012HQ Wireless</p> | <p>DR-Protectionbox DR-Skylinebox DR-Combibox LM DR-Combibox LM / AP DR-Combibox SET DR-Podo-Protectionbox DR-Cowboybox DR-Cowboy Package DR-Sleeve</p> | <p>DRKM 031012HQ DRKM 051012HQ DRKM 061012HQ 1 DRKM 061012HQ 2 DRKM 061012HQ S DRKM 071012HQ DRKM 081012HQ DRKM 081012HQ P DRKM 101012HQ</p> |

| | | | |
|---|--|---|--|
| <p>Varian</p>  | <p>Paxscan DR 4336R Leonardo DR 4336Rvet</p> | <p>DR-Protectionbox DR-Skylinebox DR-Combibox LM DR-Combibox LM / AP DR-Combibox SET DR-Podo-Protectionbox DR-Cowboybox DR-Cowboy Package DR-Sleeve</p> | <p>DRLE 03433617 DRLE 05433617 DRLE 06433617 1 DRLE 06433617 2 DRLE 06433617 S DRLE 07433617 DRLE 08433617 DRLE 08433617 P DRLE 10433617</p> |
| <p>Agfa</p>  | <p>Agfa DX-D 20</p> | <p>DR-Protectionbox DR-Skylinebox DR-Combibox LM DR-Combibox LM / AP DR-Combibox SET DR-Podo-Protectionbox DR-Cowboybox DR-Cowboy Package DR-Sleeve</p> | <p>DRAG 03DXD20 DRAG 05DXD20 DRAG 06DXD20 1 DRAG 06DXD20 2 DRAG 06DXD20 S DRAG 07DXD20 DRAG 08DXD20 DRAG 08DXD20 P DRAG 10DXD20</p> |
| <p>Agfa</p>  | <p>Agfa DX-D 25C</p> | <p>DR-Protectionbox DR-Skylinebox DR-Combibox LM DR-Combibox LM / AP DR-Combibox SET DR-Podo-Protectionbox DR-Cowboybox DR-Cowboy Package DR-Sleeve</p> | <p>DRAG 03DXD25 DRAG 05DXD25 DRAG 06DXD25 1 DRAG 06DXD25 2 DRAG 06DXD25 S DRAG 07DXD25 DRAG 08DXD25 DRAG 08DXD25 P DRAG 10DXD25</p> |
| <p>Agfa</p>  | <p>Agfa DX-D 30C Wireless</p> | <p>DR-Protectionbox DR-Skylinebox DR-Combibox LM DR-Combibox LM / AP DR-Combibox SET DR-Podo-Protectionbox DR-Cowboybox DR-Cowboy Package DR-Sleeve</p> | <p>DRAG 03DXD30 DRAG 05DXD30 DRAG 06DXD30 1 DRAG 06DXD30 2 DRAG 06DXD30 S DRAG 07DXD30 DRAG 08DXD30 DRAG 08DXD30 P DRAG 10DXD30</p> |
| <p>Agfa</p>  | <p>Agfa DX-D 35C Wireless</p> | <p>DR-Protectionbox DR-Skylinebox DR-Combibox LM DR-Combibox LM / AP DR-Combibox SET DR-Podo-Protectionbox DR-Cowboybox DR-Cowboy Package DR-Sleeve</p> | <p>DRAG 03DXD35 DRAG 05DXD35 DRAG 06DXD35 1 DRAG 06DXD35 2 DRAG 06DXD35 S DRAG 07DXD35 DRAG 08DXD35 DRAG 08DXD35 P DRAG 10DXD35</p> |
| <p>Toshiba</p>  | <p>Toshiba FDX3543RP</p> | <p>DR-Protectionbox DR-Skylinebox DR-Combibox LM DR-Combibox LM / AP DR-Combibox SET DR-Podo-Protectionbox DR-Cowboybox DR-Cowboy Package DR-Sleeve</p> | <p>DRTO 033543RP DRTO 053543RP DRTO 063543RP 1 DRTO 063543RP 2 DRTO 063543RP S DRTO 073543RP DRTO 083543RP DRTO 083543RP P DRTO 103543RP</p> |

| | | | |
|---|---|---|---|
|  | <p>Scandroid.</p> <p>Multi-functional scan stool. Ideal for the mobile veterinarian.</p> <p>Transport position.</p> <p>Position 1: 23 cm. Scanning the foot and lower parts of the leg</p> <p>Position 2: 38 cm. Scanning the front knee and middle parts of the leg</p> <p>Position 3: 51 cm. Scanning the stifle and higher parts of the leg</p> |  <p>The ideal height adjustable stool to transport to the job site or for clinic use.</p> <p>Transport measurements: ø40 x 23 cm.</p> <p>Total weight 9.3 Kg / 20.5 Lb</p> <p>The Scandroid has great stability with in total of 5 castor wheels, 2 wheels comes with locking.</p> <p>Under the seat, safely also during transport, you have some storage space for all necessary scanning tools.</p> <p>If you require we can supply the top surface of the seat with black skai leather.</p> <p>If an extra height is required we have an optional Extension with a height adjustment from 63 up to 75cm. The Extension also fits under the seat.</p>  | <p>SK 1039270</p> <p>SK 10EXT307</p> |
|  | <p>Scan Wedge.</p> <p>Simple wedge for extra extension to be placed under the hoof.</p> | <p>Material: Rubber 95° Shore</p> <p>Dimensions: 180x150x65-0mm 20° Angle</p> | <p>SW 1050065-0</p> |
|  | <p>Standoff PAD Soft</p> <p>A new soft standoff adaptable for wide range of linear-array transducers.</p> | <p>Standoff 10mm – 100 x 100 mm Standoff 20mm – 100 x 100 mm Standoff 40mm – 100 x 100 mm</p> | <p>LAS 10651010 LAS 10651030 LAS 10651040</p> |

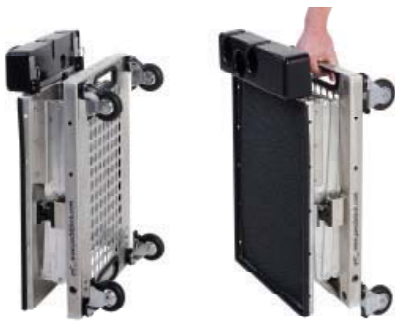

The "SCANster".

Transport position.



The optimum ultrasound solution to create the ideal working height at all times.

Lightweight durable aluminium stand.

Can be stored in vertical position to save space. Completely folded ideal for transport.



Position 1:
Height 25 cm.
For imaging the foot and pastern.



Position 2:
Height 42 cm.
For imaging the forelimbs and hind limb.



Position 3:
Height 77 cm.
For imaging the stifle, shoulder and elbow.



Position 4:
Height 112 cm.
For imaging the hip, sacrum, back, neck and endo-rectal images.

Transport: 58x41x26 cm
Total weight: 12.5 Kg
Plateau size: 55x37.5 cm
Max.load: 15 Kg

4 locking Ø80mm castor wheels.

Especially designed for the portable ultrasound.

Also works great as a clinic cart to put on small equipment to carry to the job site.

Optional with professional probe rack to store in total 2 ultrasound probes, gel bottle and [Standoff PAD](#).



VS 1045321910

CUPH 102P1B1S



The SCANster is a Class 1 European Medical Devices

MD 93/42/EEG
NL-CA002-2012-24232



The "DRagster".

Transport position.



Durable aluminium in height adjustable cart for the more heavy machines and equipment.

Standard delivered with 300N gas spring for step less variable height adjustment.

Can be stored in vertical position to save space. Easy transport position with the pull out handle.

Transport: 68x51x27 cm
Total weight: 15 Kg
Plateau size: 60 x 46 cm
Max.load: 30 Kg

2 big Ø200mm wheels
2 locking Ø125mm castors

Works great as a mobile cart for small and big equipment. So also for your DR flat panel suitcase.

VS 10322010



Position 1:
Height 25 cm.



Position 2:
Secured height 55 cm



Position 3:
Height 90 cm.



The DRagster is a
Class 1 European Medical
Devices

MD 93/42/EEG
NL-CA002-2012-24232



jiXtr.



Specifications.
 Transport:
 16 x 13 x 68 cm
 Low position:
 27 x 29 x 64 cm
 High position:
 27 x 29 x 189 cm
 Stroke range:
 150 cm
 Weight:
 4.5 Kg / 9.9 Lb

The rotation of the X-ray source in combination with the unlimited tilt of the jiXtr secures an accurate position for Skyline views.

The brand new jiXtr is yet another big step in the evolution of our products.

With the jiXtr your X-Ray source becomes an easy-to-handle, versatile and highly portable instrument that will cover every angle, every spot and every height you desire.

This light-weight, aluminium support stand allows you to work quickly and safely, and spares your back in the process.

The jiXtr can be adjusted without pre-defined increments and has a variable height-adjustment. And as for the icing on the cake: The jiXtr is compatible with all X-Ray devices that are on the market!

A "fix & release" system is integrated into the handle for easy height adjustment.

JSM 09JIXTR-14





Stat-X Vaquero.

Foldable Mobile Stand for X-ray machine.

Total weight of the Stat-X only 22 Kg.

3 Swing arm adapter positions with Quick release.

Swing arm available in 580 and 750 mm and are adjustable in height with 1 fingertip

Tailor-made bracket for X-ray unit.

Swivel range arm 5° left, 5° right.

Weight range for portable X-ray units Stat-X: 5 - 25 Kg.

Completely folded / unfolded in only 10 seconds.

External dimensions completely folded 147 x 30 x 35 cm.

Standard colour: PMS260C Satin.



Swing arm length 750 mm.
Vertical travel arm 700 mm.
Vertical travel stand 700 mm.
Max.height bracket: 207 cm.

Weight portable X-ray source:
5.0 – 15.0 Kg
15.1 – 17.0 Kg

Swing arm length 580 mm.
Vertical travel arm 500 mm.
Vertical travel stand 700 mm.
Max.height bracket: 198 cm.

Weight portable X-ray source:
17.1 – 25.0 Kg

Available for stationary use with short legs and extra contra weight to create more working space between Horse and Stat-X.

RSM 09750 LW
RSM 09750 HW

RSM 09580HW

RSM 09SLEW



The new, faster Quick-release for positioning the Swingarm in all heights



All Stat-X models are CE marked Class 1 medical product.

MD (89/392/EEG)
MDD (93/42/EEG)

Note: by ordering always provide the exact details (brand, type & weight) from your X-ray source.



Stat-X Pequeño.

Small foldable mobile stand for X-ray machine. The Stat-X Pequeño is the shorter and lighter version of the Vaquero. That way more easily to transport to the job site.

Total weight of the Pequeño only 20 Kg.

3 Swing arm adapter positions with new, faster Quick-release.

The 750mm Swing arm is adjustable in height with 1 fingertip.

Tailor-made bracket for X-ray unit.

Swivel range arm 5° left, 5° right.

Weight range for portable X-ray units.
Pequeño: 5 - 17 Kg.

Completely folded / unfolded in less than 10 seconds.

External dimensions completely folded 105 x 30 x 35 cm.

Standard colour: PMS260C Satin.



The new, faster Quick-release for positioning the Swingarm in all heights



Swing arm length 750 mm.
Vertical travel arm 700 mm.
Vertical travel stand 410 mm.
Max.height bracket: 165 cm.

Weight portable X-ray source:
5.0 – 15.0 Kg
15.1 – 17.0 Kg

PSM 09750 LW
PSM 09750 HW



All Stat-X models are CE marked Class 1 medical product.

MD (89/392/EEG)
MDD (93/42/EEG)



Stat-X
Vaquero
1.46m



Stat-X
Pequeno
1.05m

Mkr podoblock

Note: by ordering always provide the exact details (brand, type & weight) from your X-ray source.



Stat-X Espléndido.

Small foldable mobile stand for X-ray machine. The new 3D adjustable head creates a lot of extra possibilities. You can secure your generator in fixed 12° positions with a maximum of 35° left and right. When turned totally upside down the maximum height is sufficient for every X-ray exposure in your range.

Total weight of the Espléndido is only 22 Kg.

3 Swing arm adapter positions with new, faster Quick-release.

The 750mm Swing arm is adjustable in height with 1 fingertip.

New universal bracket for every X-ray unit.

Swivel range arm 5° left, 5° right.

Weight range for portable X-ray units. Espléndido: 5 - 17 Kg.

Completely folded / unfolded in less than 10 seconds.

External dimensions completely folded 114x30x35 cm.

Standard colour: PMS260C Satin.



Swing arm length 750 mm.
Vertical travel arm 700 mm.
Vertical travel stand 410 mm.
Max.height bracket: 188 cm.

Weight portable X-ray source:
5.0 – 15.0 Kg
15.1 – 17.0 Kg

The 180° rotation with fixed positions of the X-ray source secures an accurate position for Skyline views.



All Stat-X models are CE marked Class 1 medical product.

MD (89/392/EEG)
MDD (93/42/EEG)

ESM 09750 LW
ESM 09750 HW

Note: by ordering always provide the exact details (brand, type & weight) from your X-ray source.



Stat-X Bucky.

Foldable Mobile Stand with 3D adjustable and mountable cassette holder.

Total Weight of the Stat-X only 22 Kg.

3 Swing arm adapter positions with Quick release.

The 750mm Swing arm is adjustable in height with only 1 fingertip.

Cassette holder is 3D adjustable.

Inclusive tailor-made Cassette holder.

Swivel range arm 5° left, 5° right.

Weight range for Bucky: 3 - 15 Kg.

Completely folded / unfolded in only 10 seconds.

Standard colour: PMS260C Satin.



The new, faster Quick-release for positioning the Swingarm in all heights.



Mobile Bucky inclusive a Cassette holder for Analogue, CR or DR system. Optional tailor-made grid compartment.

Swing arm length 750 mm.
Vertical travel arm 700 mm.
Vertical travel stand 700 mm.

Weight X-ray Cassette / Panel:
3.0 – 15.0 Kg

Available for stationary use with short legs and extra contra weight to create more working space between Horse and Stat-X.

Optional digital protractor: Aiangle



All Stat-X models are CE marked Class 1 medical product.

MD (89/392/EEG)
MDD (93/42/EEG)

BSM 09750 LW

RSM 09SLEW

AI 10RAS2C

Note: by ordering always provide the exact details (brand, measurements & weight) from your used system.



Stat-X Infinity.

The Stat-X Infinity is the ultimate foldable stand.

It has all the options needed to x-ray whatever part needed quickly and safely.

Total weight of the Stat-X only 26 Kg.

3 Swing arm adapter positions with new Quick release.

Swing arm available in 580 and 750 mm and are adjustable in height with 1 fingertip

Tailor made bracket for X-ray unit.

Weight range for portable X-ray units Stat-X: 5 - 25 Kg.

Weight range for Cassette holder Infinity: 3 - 17 Kg.

Cassette holder is 3D adjustable.

Tailor-made Cassette holder.

Swivel range arm 5° left, 5° right.

Completely folded / unfolded in only 10 seconds.

Standard colour: PMS260C Satin.

Provided with the new, faster Quick-release for positioning the Swingarm in all heights.



Swing arm length 750 mm.
Vertical travel arm 700 mm.
Vertical travel stand 700 mm.
Max.height bracket: 207 cm.

Weight portable X-ray source:

5.0 – 15.0 Kg
15.1 – 17.0 Kg

Swing arm length 580 mm.
Vertical travel arm 500 mm.
Vertical travel stand 700 mm.
Max.height bracket: 198 cm.

Weight portable X-ray source:
17.1 – 25.0 Kg

Available for stationary use with short legs and extra contra weight to create more working space between Horse and Stat-X.

Optional digital protractor: Aiangle.



All Stat-X models are CE marked Class 1 medical product.

MD (89/392/EEG)
MDD (93/42/EEG)

ISM 09750 LW
ISM 09750 HW






ISM 09580HW

RSM 09SLEW

AI 10RAS2C

Note: by ordering always provide the exact details (brand, type & weight) from your complete X-ray system.

| | | | |
|---|--|---|---|
|     | <p>X-Ray clip. Stainless steel clip for analogue and CR-Cassettes.</p> <p>Sarah marker. Marker with integrated libel to be always sure about the horizontal plane.</p> <p>Multimarker. Stainless steel all in one marker with suction cup.</p> <p>Buckymarker. Marker with suction cup for your cassette holder to provide information about the horizontal plane.</p> | <p>L / R LF / RF LV / RV LH / RH LAT / MED</p>  <p>L / R with suction cup</p> <p>LF / RF / LH / RH LV / RV / LH / RH LAT / MED / DV / VD</p> <p>Stainless steel. 360° rotation. 48 x 44 x 4 mm</p> | <p>RC 10LR18 RC 10LFRF18 RC 10LVVRV18 RC 10LHRH18 RC 10LATMED18</p> <p>SC 10LRZN</p> <p>RCM 10ENZN RCM 10DEZN RCM 10LMZN</p> <p>MKB 1048444</p> |
|  | <p>Warning sign. Mobile warning sign that radiation will be released.</p> | <p>Language Dutch Language English Language German Language French</p> | <p>WS 10600NL WS 10600EN WS 10600DE WS 10600FR</p> |
|  | <p>Hoofplaster. For shadow free hoof recordings.</p> | <p>5 x 250 gram 25 x 250 gram 50 x 250 gram</p> | <p>HP 10BL25005 HP 10BL25025 HP 10BL25050</p> |
|  | <p>VET Tools. To work quick and safe.</p> | <p>Apron "VET" Adjustable Hooftester</p> <p>VET Tools Set: Apron "VET" Soft faced hamer Clench cutter Straight blade Loop Knife Nail puller Pull off Adjustable Hooftester Canvas Storage bag</p> | <p>VV 10VET VV 10WERK07</p> <p>VV 10SET</p> |
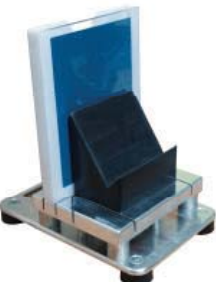
| | | | |
|--|--|---|--|
|  | <p>Aluminium Wedge.</p> <p>Filter the radiation with high purity homogeneous quality aluminium for better back / neck X-rays.</p> | <p>Small wedge with Velcro for use on your X-ray generator. Dimensions: 195x75mm / 23°</p> <p>Small wedge with Velcro for use on your X-ray generator. Dimensions: 195x50 mm / 45°</p> <p>Large concave wedge by Prof. J.M. Denoix. Provide with a contra weight for easy fixation in every position. Dimensions: 400x200x40-0 mm</p> | <p>AK 101957530</p> <p>AK 101955045</p> <p>AK 1040025040</p> |
|  | <p>Aiangle.</p> <p>Great device to measure the exact right angle. Easy to place on your Generator and cassette holder.</p> | <p>Measuring range: +/- 180° Resolution: 0.1° Magnetic base. Tilt display 0°-90° Included 2AAA batteries.</p> | <p>AI 10RAS2C</p> |
|  | <p>Grid compartment.</p> | <p>Tailor-made Grid compartment for Podoholder and DR flat panel products.</p> | <p>MR 10RO01</p> |
| <div data-bbox="140 1288 406 1489">  </div> <p>Podoblock + Podoholder with Grid compartment.</p> <div data-bbox="167 1724 383 2004">  </div> | <p>Aluminium Grid.</p> <p>Aluminium cover. Parallel lines. 70 lines / cm. 178 lines / inch. Focus: Infinity. Ratio 8:1.</p> <p>Carbon Grid.</p> <p>Carbon cover. Fibre interspaces. Parallel lines 70 lines / cm. 178 lines / inch. Focus: Infinity. Ratio 13:1.</p> | <p>18x24 Film.size: 180x240mm 18x24 Cass.size: 200x260mm 24x30 Film.size: 240x300mm 24x30 Cass.size: 260x320mm 35x43 Film size: 356x432mm</p> <p>18x24 Film.size: 180x240mm 18x24 Cass.size: 200x260mm 24x30 Film.size: 240x300mm 24x30 Cass.size: 260x320mm 35x43 Film size: 350x430mm</p> | <p>RL 10182415 RL 10202615 RL 10243015 RL 10263215 RL 10354315</p> <p>RLC 10182415 RLC 10202615 RLC 10243015 RLC 10263215 RLC 10354315</p> |

All other Grid specifications (Focussed, Lines/cm, Ratio) are available on request.

| | | | |
|---|---|---|--|
|  | <p>X-ray Apron.</p> <p>Lightweight nylon. 0.50 mm Pb Front 0.25 mm Pb Back. Snap fasteners.</p>  | <p>Small 50 x 110 cm Medium 60 x 120 cm Large. 75 x 120 cm</p> <p>For available colours and prints see table next page.</p> | <p>RJ 10CLI05S RJ 10CLI05M RJ 10CLI05L</p> |
|  | <p>X-ray Apron “de Luxe”</p> <p>Lightweight nylon with integrated Thyroid shield protection. 0.50mm Pb Front. 0.25mm Pb Back. 0.50mm Pb Thyroid. Snap fasteners.</p>  | <p>Small 50 x 110 cm Medium 60 x 120 cm Large. 75 x 120 cm</p> <p>For available colours and prints see table next page.</p> | <p>RJ 10CLI05STS RJ 10CLI05MTS RJ 10CLI05LTS</p> |
|  | <p>X-ray glove.</p> <p>Flexible vinyl glove with 0.5mm Pb.</p> | <p>Pb 0.5 mm Size Small Pb 0.5 mm Size Medium Pb 0.5 mm Size Large</p> | <p>RH 10LF05S RH 10LF05M RH 10LF05L</p> |
|  | <p>Storage.</p> <p>Apron rack Glove rack</p> | <p>5 Garments, wall mount. 3 Pair gloves, wall mount.</p> | <p>RR 10CH05 RR 10GL03</p> |
|  | <p>Thyroid shield.</p> <p>Lightweight Nylon. Velcro fastener.</p> | <p>One size fits all. Pb 0.5mm Colour: Blue</p> | <p>RS 10BL05M</p> |
|  | <p>X-ray Protection Goggles.</p> | <p>Pb 0.5mm</p> | <p>RB 10OP05</p> |

Measurements X-ray cassettes. Analogue & Digital CR.

Please check the Cassette measurements if you are not sure what type of cassette you have. We supply cassette holders for every X-ray system. If you have different measurements cassettes, please contact us.

| Film size (cm) | Cassette size (mm) | Film size (inch) | MXR Article number |
|----------------|--------------------|------------------|--------------------|
| 18 x 24 | 208 x 267 | |1824... |
| 24 x 30 | 267 x 327 | |2430... |
| 35 x 35 | 385 x 385 | 14" x 14" |3535... |
| 35 x 43 | 385 x 460 | 14" x 17" |3543... |
| 43 x 43 | 460 x 460 | 17" x 17" |4343... |
| 20 x 25 | 230 x 280 | 8" x 10" |2328... |
| 25 x 30 | 280 x 330 | 10" x 12" |2833... |
| 25 x 35 | 280 x 390 | 10" x 14" |1014... |
| 20 x 30 | 230 x 330 | 8" x 12" |2030... |
| 20 x 40 | 230 x 425 | 8" x 16 " |2040... |
| 30 x 40 | 330 x 425 | 12" x 16 " |3040... |

Colour card X-ray Apron and X-ray Apron "de Luxe"

All colours available without extra cost.

Place the requested colour number on your order behind the Apron article number.




Nr. 8 Dark Gray



Nr. 2-1 Red with Horse




Nr. 13 Dark Blue

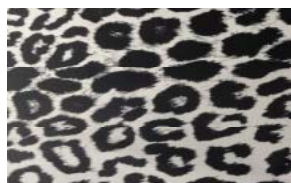

Nr. 12 Leopard




Nr. 6-1 Dark Red

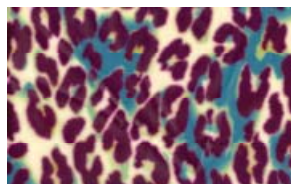

Nr. 12-3 Purple Leopard



Nr. 1 Dark Green

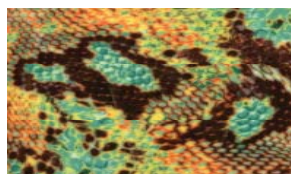

Nr. 5-5 Python



Nr. 4 Dark Purple



Nr. 4-2 Purple with Motif



Blue Cape Vet.

This apron is especially developed for those veterinarians that work with a multitude of animals, either inside or outside.

It is a pleasure to wear under all circumstances and extremely easy to clean.

Consider these 10 advantages of this state-of-the-art apron and don't waste your time shopping around for a better deal.

1. Elastic waist band
2. Detachable sleeves
3. Semi-elastic cuffs
4. Backside skirt system
5. Two outer and inner zippable pockets
6. Wind tight
7. Water repellent
8. Wrinkle-free
9. Light-weight
10. Comfortable to wear

Female - XS
 Female - S
 Female - M
 Female - L
 Female - XL
 Female - XXL

BCP 10FEMALE-XS
 BCP 10FEMALE-S
 BCP 10FEMALE-M
 BCP 10FEMALE-L
 BCP 10FEMALE-XL
 BCP 10FEMALE-XXL



Male - XS
 Male - S
 Male - M
 Male - L
 Male - XL
 Male - XXL

BCP 10MALE-XS
 BCP 10MALE-S
 BCP 10MALE-M
 BCP 10MALE-L
 BCP 10MALE-XL
 BCP 10MALE-XXL



| | | |
|---|--|--------------------|
|  | <p>Digital Hoof Tester.</p> <p>The Digital Hoof Tester™ measures the sensitivity of a horse's hoof. The main difference is that this is the first digital hoof tester. It is made from lightweight aluminum and has stainless steel side reinforced plates.</p> <p>Mechanical features:</p> <ul style="list-style-type: none"> * Hoof contact swivel pads * Center axis with return spring * Modular, ergonomic design * Smooth contours * Single-hand grip <p>Digital Features:</p> <ul style="list-style-type: none"> * Two digit LCD readout * Set point with high/low reply * Automatic shut off * Invertible display * Low power / long battery life | <p>VV 10DIG07</p> |
|  | <p>Digital Hoof Gauge.</p> <p>The Digital Hoof Gauge™ accurately measures the angle of a horse's hoof up to 0.1th of a degree.</p> <p>Collapsible, Rugged, Time saving, Auto On/Off, Water tight, Accurate, Long lasting battery</p> | <p>DHG 10PEP01</p> |
|  | <p>123Shoeing Kit.</p> <ol style="list-style-type: none"> 1. Shoe Puller 2. Hammer 3. Clinch Cutter <p>The 123Shoeing Kit is an essential tool for de-shoeing a horse. Whether you're an equine Vet or a Farrier. You should take this high quality, professional 3-in-1-tool in your basic kit. Use this quality farrier-hammer with the clinch cutter for removing nails during de-shoeing. You can use this indispensable tool to cut nail clinches when removing shoes and to remove "hidden" nails embedded in hoof wall, which can probably ruin your nipper.</p> <p>The essential tool for practical minded professionals!</p> | <p>PZT 10123SK</p> |



MXR Podoblock B.V.

Esweg 31
9475 TM
Midlaren
The Netherlands

| | |
|--------------|------------------|
| Nico Gorter | 0031-653 794 282 |
| Alex Berends | 0031-651 828 117 |
| Office | 0031-508 200 257 |
| Fax | 0031-508 200 212 |

| | |
|----|--|
| I: | www.podoblock.com |
| E: | info@mxr.nl |
| E: | sales@mxr.nl |