



CMI

Calico Medical


INCORPORATED

Zenith® Proudly Supports Chiropractic Worldwide For Nearly A Century




Zenith® 230-Plus

Zenith® 225 Diversified Hylo

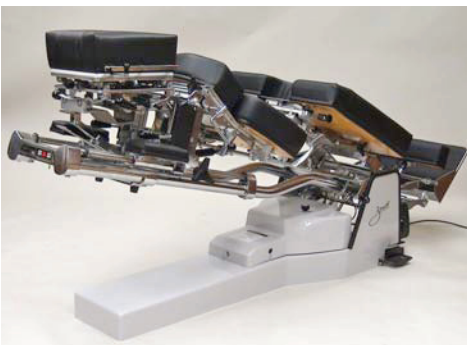


Zenith® E44 THOMPSON®

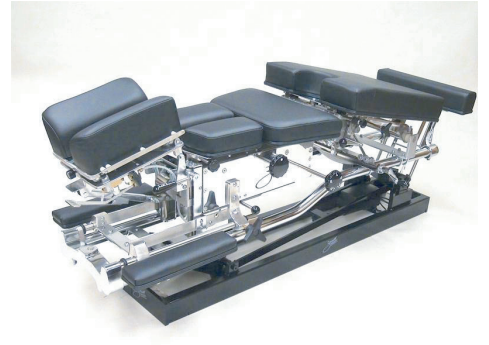

Zenith®

Zenith® 125 Toggle

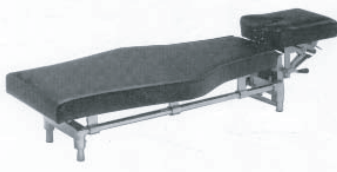

1789 Notre Dame Dr.
St Agatha ON. N0B 2L0
Phone: 519.885.3254
Email: office@calicomedical.com
www.calicomedical.com

Zenith® 100 F/D





CMI
Calico Medical
 INCORPORATED

Zenith® FAMILY OF ADJUSTING TABLES

Developed by Clay Thompson, DC., *Zenith* Thompson Techniques used by Thompson® adjusting table with pneumatically driven, segmented drop system, which quickly lowers the selection of the patient's body corresponding with the spinal region being adjusted. The doctor thrusts at high speed, using minimal force, because while the thrust initiates movement, the fast drop carries the joint through the remainder of its range of motion.

With more chiropractic tables operating in more colleges and universities than any other single supplier, or even with all other suppliers combined, the *Zenith*® family of adjusting tables is built to exacting specifications to arrive at a solution that is rugged, dependable and aesthetically beautiful.

Zenith® quality is backed by our **EXCLUSIVE 5 YEAR LIMITED PARTS WARRANTY**, that no other company has ever tried to match.



Zenith® STATIONARY

Rugged, dependable, affordable and available in no drop through 4-drop *Zenith*® THOMPSON® models.



Zenith® ELEVATION-SERIES™

For the doctor that chooses an elevation table, the *Zenith*® Elevation-Series provides the total solution.



Zenith® HYLO & VERTI-LIFT

The most popular chiropractic adjusting solution in the world. The right choice for comfortable patient positioning.